**PORTUGUÊS**

Implante

**INGLÊS**

Implant

**ESPAÑHOL**

Implante

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Technical name: Dental Implants (Osseointegrable)  
Trade name: implant

**ATTENTION!!!**

**SINGLE USE PRODUCT. REPROCESSING PROHIBITED.**

**DO NOT USE IF PACKAGING IS DAMAGED OR BROKEN.**

**STERILE PRODUCT. DO NOT RESTERILIZE. STERILIZED BY GAMMA RADIATION.**

**MANUFACTURER RECOMMENDS SINGLE USE.**

**Business models:**

Implant Stretto  
Implant ALTIVUS CM20  
Implant ALTIVUS Z CM20  
Implant ALTIVUS Z TL CM20  
Implant Veloce CM  
Implant Veloce HE  
Implant Veloce HI  
Implant Conical HI  
Implant Conical HE  
Implant Cone Morse  
Implant Collo HE  
Implant Conical HE IL Maestro  
Implant Conical HI IL Maestro  
Implant Due Cone IL Maestro  
Implant Cylindrical HI  
Implant Cylindrical HE

Note: To consult the product reference/code, see Annex I.

**1- FORMS OF COMMERCIAL PRESENTATION OF THE PRODUCT**

Implacil Implants are presented in individual packaging. The implant is packaged inside a capsule that is inserted into a blister sealed with Tyvek. The blister is placed inside a cardboard box.

List of items contained inside the package:

Models: Implant Cylindrical HI, Implant Cylindrical HE.  
Implant Conical HI, Implant Conical HE, Implant Cone Morse  
- 01 Implant  
- 01 Implant Cover  
- 01 Assembler  
- 01 Assembler Screw

Models: Implant Stretto, Implant Veloce CM, Implant Veloce HE, Implant Veloce HI, Implant Conical HE IL Maestro, Implant Conical HI IL Maestro, Implant Due Cone IL Maestro, Implant Collo HE, Implant ALTIVUS CM20, Implant ALTIVUS Z CM20 and Implant ALTIVUS Z TL CM20  
- 01 Implant  
- 01 Implant Cover

This instruction for use will be available on the Implacil Osstem website (<https://implacilosstem.com.br/>) in

accordance with current legislation on the provision of instructions for use in non-printed format for health products – RDC 751/2022. Printed instructions for use can be requested through Customer Service – 0800.779.0011, and will be sent at no additional cost.

To locate it, check the “IFU” number and the version of the instruction that is indicated on the Label of the purchased product.

To read the instructions for use on the website, you need Adobe Acrobat Reader.

User Alert: Always observe the correlation of the version of the instructions for use indicated on the label of the purchased product, as made available by Implacil Osstem.

**2- DETAILED DESCRIPTION**

**2.1 – Commercial Model: Implant Stretto**

The IMPLANT STRETTO consists of a conical body. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring, increasing the osseointegration surface.

The implant is indicated for little prosthetic space or little bone structure and should only be used in the central and lateral incisors and canines.

On the longitudinal axis of the upper part of the implant there is a central thread (M1.4) that allows the prosthetic component to be fixed by a screw. It is produced by turning titanium bars in accordance with ASTM F-67 standard. It is offered in diameters 3.0 and in lengths 10, 11.5, 13 and 15mm.

**2.2- Commercial Model: Implante ALTIVUS CM20**

The ALTIVUS IMPLANT CM20 consists of a cylindrical body. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchorage by increasing the osseointegration surface.

The implant is indicated for rehabilitation in patients with atrophic jaws, offering a safe and efficient alternative when bone grafts or the installation of shorter implants are not

feasible. Its design is optimized to provide maximum primary stability, bone anchorage in the transnasal or pterygoid region, and ease of surgical protocol, and is indicated in conjunction with CM components.

On the longitudinal axis of the upper part of the implant, there is a central thread (M1.8) that allows the prosthetic component to be fixed with a screw.

It is offered in diameters 3.75 and 4.0mm; and in lengths 19; 21; 23 and 25mm.

### **2.3- Commercial Model: Implante ALTIVUS Z CM20**

The Altivus Implant Z CM20 (Fully Sandblasted Neck) has a cylindrical body made of titanium. The entire surface of the implant, including the neck region, is mechanically sandblasted, expanding the contact area and promoting anchorage and osseointegration along its entire length.

This implant is indicated for the rehabilitation of patients with atrophic jaws, representing a safe and effective alternative in cases where bone grafts or shorter implants are not feasible. Its design has been optimized to provide maximum primary stability, promote anchorage in the zygomatic region, and simplify the surgical protocol. Furthermore, its geometry contributes to increased implant resistance in areas of low bone density. It is recommended for use in conjunction with the CM and Z Abutment components.

The upper portion of the implant features a central thread (M1.8) along the longitudinal axis, allowing fixation of the prosthetic component with a screw.

It is offered in diameters 3.75 and 4.0mm; and in lengths 19; 21; 23 and 25mm.

### **2.4- Commercial Model: Implante ALTIVUS Z TL CM20**

The Implant Altivus Z TL CM20 has a cylindrical body made of titanium. Its surface is mechanically treated by sandblasting, which expands the contact area, favors anchorage, and enhances osseointegration. Only the neck region maintains a low-roughness finish, specially designed to optimize adaptation to bone tissue.

This implant is indicated for the rehabilitation of patients with atrophic jaws, representing a safe and effective alternative in cases where bone grafts or the installation of shorter implants are not feasible. Its design has been optimized to provide maximum primary stability, promote anchorage in the zygomatic region, and simplify the surgical protocol. Furthermore, its geometry contributes to increased implant resistance in areas of low bone density. It is recommended for use in conjunction with CM and Z-Abutment components.

At the top, the implant features a central thread (M1.8) along the longitudinal axis, allowing fixation of the prosthetic component with a screw.

It is offered in diameters 3.75 and 4.0mm; and in lengths 19/21/23 and 25mm.

### **2.5 – Commercial Model: Implant Veloce CM**

The IMPLANT VELOCE CM consists of a cylindrical body with double threads that aid its installation and anchoring in the bone bed along its entire length. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring, increasing the osseointegration surface. The double threads serve to facilitate and speed up its insertion.

On the longitudinal axis of the upper part of the implant there is a central thread that allows the prosthetic component to be fixed by a screw. The design of this implant favors load distribution, and is mainly indicated for use in low-density bone. It is produced by turning titanium bars in accordance with ASTM F-67 standard.

It is offered in diameters 3.3; 3.75; 4.3 and 4.75mm; and in lengths 8; 10; 11.5; 13; 15 and 17mm.

### **2.6 – Commercial Model: Implant Veloce HE**

The IMPLANT VELOCE HE consists of a cylindrical body with double threads along its entire length that aid its installation and anchoring in the bone bed. The implant body is mechanically treated (sandblasting), a process that

increases its contact surface and facilitates anchorage increasing the osseointegration surface. The double threads serve to facilitate and speed up insertion.

The external hexagonal connection is a well-established prosthetic fitting that has been used worldwide for over 40 years. There is a central thread on the longitudinal axis of the upper part of the implant that allows the prosthetic component to be fixed with a screw. The design of this implant favors load distribution, and is mainly indicated for use in low-density bone or even for immediate loading. It is produced by turning titanium bars in accordance with ASTM F-67.

It is offered in diameters 3.3; 3.75; 4.3 and 4.75mm; and in lengths 8; 10; 11.5; 13; 15 and 17mm.

### **2.7 – Commercial Model: Implant Veloce HI**

The IMPLANT VELOCE HI consists of a cylindrical body with double threads along its entire length that aid its installation and anchoring in the bone bed. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring, increasing the osseointegration surface. The double threads serve to facilitate and speed up its insertion.

The internal hexagonal connection makes this system one of the simplest to handle. On the longitudinal axis of the upper part of the implant there is a central thread (g) that allows the prosthetic component to be fixed by a screw. The design of this implant favors load distribution, being mainly indicated for use in low-density bone or even for immediate loading. It is produced by turning titanium bars according to ASTM F-67 standard.

It is offered in diameters 3.3; 3.75; 4.3 and 4.75mm; and in lengths 8; 10; 11.5; 13; 15 and 17mm.

### **2.8 – Commercial Model: Implant Conical HI**

The IMPLANT HI Conical consists of a conical body with threads along its entire length that aid its installation and anchoring in the bone bed. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring by

increasing the osseointegration surface, maintaining only in the neck region a low-roughness titanium finish to improve its adhesion to the bone tissue. In the upper part, this implant has micro-spirals, which improve distribution in the bone cortex. The implant also has three millings at the apex, distributed evenly (120°) in order to facilitate its insertion. It is produced with a slightly divergent neck. The internal hexagonal connection makes this system one of the simplest to handle. In the longitudinal axis of the upper part of the implant there is a central thread that allows the prosthetic component to be fixed by a screw. The design of this implant favors load distribution, and is mainly indicated for use in low-density bone or even for immediate loading. It is produced by turning commercially pure titanium bars, meeting ASTM F-67 and ASTM F136 standards.

It is offered in diameters of 3.50mm; 4.00mm and 5.0mm; and in lengths of 7mm, 8mm, 9mm, 10mm, 11mm, 12mm, 13mm, 14mm, 15mm, 16mm, 17mm and 18mm.

### **2.9 – Commercial Model: Implant Conical HE**

The IMPLANT Conical HE– External Hexagon – consists of a conical body with threads along its entire length that aid its installation and anchoring in the bone bed. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring by increasing the osseointegration surface, maintaining only in the neck region a low-roughness finish of the titanium to improve its adhesion to the bone tissue. Also in the upper part, this implant has micro-spirals, which improve distribution in the bone cortex. The implant also has three millings at the apex, distributed evenly (120°) in order to facilitate its insertion. It is produced with a slightly divergent neck. The external hexagonal connection is a well-established prosthetic fitting, used worldwide for over 40 years. On the longitudinal axis of the upper part of the implant there is a central thread that allows the prosthetic component to be fixed by a screw. The design of this implant favors load distribution, and is mainly indicated for use in low-density bone or even for immediate loading. It is produced by turning commercially pure titanium bars, meeting the ASTM F-67 and ASTM F136 standards. It is offered in diameters of 3.50 mm; 4.00 mm and 5.0 mm; and in lengths of 7 mm, 8 mm, 9 mm, 10 mm, 11 mm, 12 mm, 13 mm, 14 mm, 15 mm, 16 mm, 17 mm and 18 mm. It can be offered with an assembler. The HE Conical IMPLANT combines the efficiency of the external hexagonal prosthetic fitting, recognized worldwide, with the practicality provided by the conical design.

It is offered in diameters of 3.50mm; 4.00mm and 5.0mm; and in lengths of 7mm, 8mm, 9mm, 10mm, 11mm, 12mm, 13mm, 14mm, 15mm, 16mm, 17mm and 18mm.

#### **2.10 – Commercial Model: Implant Cone Morse**

The IMPLANT CONE MORSE consists of a conical body with threads that aid its installation and anchoring in the bone bed along its entire length. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring, increasing the osseointegration surface. This surface treatment is performed on the entire implant, since it is indicated for installation 1 to 2 mm below the bone crest. The implant also has three millings at the apex, distributed evenly (120°) to facilitate its insertion. It is produced with a slightly divergent neck. The Morse Taper connection of this implant favors its use for unitary elements. On the longitudinal axis of the upper part of the implant there is a central thread that allows the fixation of the prosthetic component by a screw. The design of this implant favors load distribution, being mainly indicated for use in low-density bone. It is produced from the turning of commercially pure titanium bars, meeting ASTM F-67; and ASTM F136 standards.

It is offered in diameters of 3.50mm; 4.00mm, 4.5mm and 5.0mm; and in lengths of 8mm, 9mm, 10mm, 11mm, 12mm, 13mm, 14mm, 15mm, 16mm, 17mm and 18mm.

#### **2.11 – Commercial Model: Implant Collo HE**

IMPLANT Collo HE consists of a cylindrical body with threads that aid its installation and anchoring in the bone bed. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring by increasing the osseointegration surface, maintaining only in the neck region a low roughness finish of the titanium to improve its adhesion to the bone tissue. The elongated neck with this smoothing in the upper profile was designed to eliminate the need for implantoplasty.


The implant has three millings at the apex, distributed evenly (120°) to facilitate its insertion. Its seating platform has a diameter of 4.0 mm.

The external hexagonal connection is a well-established prosthetic fitting that has been used worldwide for over 40 years. On the longitudinal axis of the upper part of the implant there is a central thread that allows the prosthetic component to be fixed by a screw. It is produced by turning commercially pure titanium bars, meeting the ASTM F-67 standard;

The IMPLANT Collo HE has the prosthetic fitting that has been the leading choice of professionals around the world for decades.

It is offered in a diameter of 4.0mm; in lengths of 8mm, 10mm, 11.5mm, 13mm, 15mm and 18mm.


#### **2.12 – Commercial Model: Implant Conical HE IL Maestro**

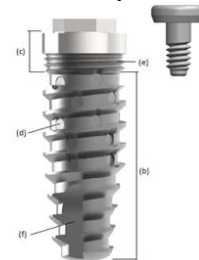The IMPLANT CONICAL HE IL MAESTRO consists of a conical body with threads along its entire length that aid its installation and anchoring in the bone bed. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring by increasing the osseointegration surface, maintaining only in the neck region a low-roughness titanium finish to improve its adhesion to the bone tissue. Still aiming at improving osseointegration, this model has radial channels that differentiate this implant profile from similar models that are already registered, as these channels provide better integration; such channels function as pockets of autogenous bone residues, reducing the implant placement torque, inducing an improvement in integration time. Still in the upper part, this implant has micro-spirals, which improve distribution in the bone cortex. The implant also has three millings at the apex, distributed evenly (120°) to facilitate its insertion. It is produced with a slightly divergent neck. The external hexagonal connection is a well-established prosthetic fitting, used worldwide for over 40 years. On the longitudinal axis of the upper part of the implant there is a central thread that allows the prosthetic component to be fixed by a screw. The design of this implant favors load distribution, being mainly indicated for use in low-density bone or even for immediate loading. It is produced by turning titanium bars according to ASTM F-67 standard.

It is offered in diameters 3.50; 4.00 and 5.0mm; and in lengths 5; 6; 7; 8; 9; 10; 11; 12; 13; 14 and 15mm.

### 2.13 – Commercial Model: Implant Conical HI IL Maestro

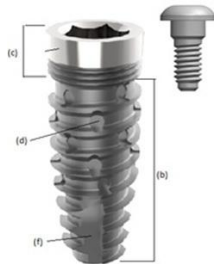The IMPLANT CONICAL HI IL MAESTRO consists of a conical body with threads along its entire length that aid its installation and anchoring in the bone bed. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring by increasing the osseointegration surface, maintaining only in the neck region a low-roughness titanium finish to improve its adhesion to the bone tissue. Also aiming at improving osseointegration, this model has radial channels that differentiate this implant profile from similar models that are already registered, as these channels provide better integration; such channels function as pockets of autogenous bone residues, reducing the implant placement torque, inducing an improvement in integration time. In the upper part, this implant has micro-spirals, which improve distribution in the bone cortex. The implant also has three millings at the apex, distributed evenly ( $120^\circ$ ) to facilitate its insertion. It is produced with a slightly divergent neck. The internal hexagonal connection makes this system one of the easiest to handle. On the longitudinal axis of the upper part of the implant there is a central thread that allows the prosthetic component to be fixed by a screw. The design of this implant favors load distribution, being mainly indicated for use in low-density bone or even for immediate loading. It is produced by turning titanium bars according to ASTM F-67 standard.

It is offered in diameters 3.50; 4.00 and 5.0mm; and in lengths 5; 6; 7; 8; 9; 10; 11; 12; 13; 14 and 15mm.

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The IMPLANT DUE CONE IL MAESTRO consists of a conical body with threads that aid its installation and anchoring in the bone bed along its entire length. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring, increasing the osseointegration surface. Also aiming at improving osseointegration, this model has radial channels that differentiate this implant profile from similar models that are already registered, as these channels provide better integration; such channels function as pockets of

autogenous bone residue, reducing the implant placement torque, inducing an improvement in integration time. This surface treatment is performed on the entire implant, since it is indicated for installation 1 to 2 mm below the bone crest. The implant also has three millings at the apex, distributed evenly ( $120^\circ$ ) in order to facilitate its insertion. It is produced with a slightly divergent neck. The Morse Taper connection of this implant favors its use for single elements. It has an internal 2.00 mm hexagon used both for insertion and as an index for prosthetic elements. On the longitudinal axis of the upper part of the implant there is a central thread that allows the prosthetic component to be fixed by a screw. The design of this implant favors load distribution, being mainly indicated for use in low-density bone. It is produced from the turning of titanium bars according to ASTM F-67 standard.

It is offered in diameters 3.50; 4.00; 4.50 and 5.0mm; and in lengths 5; 6; 7; 8; 9; 10; 11; 12; 13; 14 and 15mm.

### 2.15- Commercial Model: Implant Cylindrical HI

The Implant Cylindrical HI consists of a cylindrical body with threads that aid its installation and anchoring in the bone bed along its entire length. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring by increasing the osseointegration surface, maintaining only in the neck region a low-roughness titanium finish to improve its adhesion to the bone tissue. The implant has three millings at the apex, distributed evenly ( $120^\circ$ ) to facilitate its insertion. Its seating platform has the following diameters: 3.50 mm, which corresponds to the 3.30 mm implant; 4.00 mm, which corresponds to the 3.75 mm and 4.30 mm implant; 5.00 mm, which corresponds to the 4.75 mm implant. It is produced with a slightly divergent neck and a self-threading conical apex. The internal hexagonal connection makes this system one of the simplest to manipulate. On the longitudinal axis of the upper part of the implant there is a central thread that allows the prosthetic component to be fixed by a screw. It is produced by turning commercially pure titanium bars, meeting the ASTM F-67 and ASTM F136 standards.

It is offered in diameters 3.30mm; 3.75mm; 4.30mm and 4.75mm; and in lengths 7.3mm, 8mm, 9mm, 10mm, 11mm, 12mm, 13mm, 14mm, 15mm, 16mm, 17mm and 18mm.

## **2.16- Commercial Model: Implant Cylindrical HE**

The Implant Cylindrical HE consists of a cylindrical body with threads that aid its installation and anchoring in the bone bed along its entire length. The implant body is mechanically treated (sandblasted), a process that increases its contact surface and facilitates anchoring by increasing the osseointegration surface, maintaining only in the neck region a low-roughness titanium finish to improve its adhesion to the bone tissue. The implant has three millings at the apex, distributed evenly (120°) to facilitate its insertion. Its seating platform has the following diameters: 3.50 mm, which corresponds to the 3.30 mm implant; 4.00 mm, which corresponds to the 3.75 mm and 4.30 mm implant; 5.00 mm, which corresponds to the 4.75 mm implant. It is produced with a slightly divergent neck and a self-threading conical apex. The external hexagonal connection is a well-established prosthetic fitting, used worldwide for over 40 years. On the longitudinal axis of the upper part of the implant there is a central thread that allows the prosthetic component to be fixed by a screw. It is produced by turning commercially pure titanium bars, meeting the ASTM F-67 and ASTM F136 standards.

The HE Cylindrical IMPLANT has the prosthetic fitting that has been the leading choice of professionals worldwide for decades. Its neck, or counter-sink, adapts well to the cortical bone, efficiently distributing the occlusal load.

It is offered in diameters of 3.30mm; 3.75mm; 4.0mm and 4.75mm; and in lengths of 8mm, 10mm, 11.5mm, 13mm, 15mm and 18mm.

Note: For information on implant compatibility, see Annex II and for the sequence of burs, see Annex III.

### **3- INDICATION OF USE AND PURPOSE**

Indication: Patients who require total or partial oral rehabilitation.

Purpose: To perform the function of the root of the missing tooth.

### **4- PHYSICAL PRINCIPLE AND FUNDAMENTALS OF PRODUCT TECHNOLOGY, APPLIED TO ITS OPERATION AND ACTION**

The operating principle of a dental implant is the biocompatibility of the material from which it is made, which allows its bone integration (maxillary and mandibular), performing the function of the root of the lost tooth. The thread is designed to conveniently distribute the masticatory load. The porosity of the surface exerts bone

retention, which is distributed homogeneously in the surface area and also acts in the microdistribution of the masticatory load. The apical recesses act after osseointegration as mechanical retention.

The implant is inserted into the bone through an intraoral surgical process, fixed by locking in bone type I, II, III or IV. IMPLACIL implants were designed for single, multiple, total restorations and also for the Immediate Loading system.

### **5- HOW TO USE THE PRODUCT**

- Open the outer packaging. Check that the blister is perfectly sealed. Reject the product if the packaging is damaged.
- Start the sequence of the drills (see appendix III) to create the cavity where the implant will be inserted. Open the blister only at the time of insertion of the implant.
- Remove the capsule from the blister and open it carefully. Do not throw away the cap, as the implant cap is inside it.
- Remove the implant from inside the capsule using the key.
- Insert the implant into the newly drilled cavity. Maximum insertion speed of 20 RPM. If done manually, the maximum torque to be used is 60 N.cm for the Cylindrical HI/HE, Conical HI/HE, Cone Morse, Veloce CM/HE/HI and ALTIVUS CM20/Z CM20/ Z TL CM20 models, and 40 N.cm for the Conical HI/HE IL Maestro and Due Cone IL Maestro models.
- Once inserted, remove the mounter (when applicable) using a 1.17mm hex key.
- Remove the implant cap from the capsule lid.
- Seal with the implant cover.
- Installation complete.

### **6- COMPOSITION**

As raw material for implants, it is produced from the turning of titanium bars – ASTM F67.

### **7- VALIDITY PERIOD**

05 years from the date of manufacture.

### **8- STORAGE, TRANSPORTATION AND HANDLING CONDITIONS**

Storage: Keep in a place protected from rain and direct sunlight, and in its original packaging.

Transportation: When transporting, protect the packaging from direct rain.

Packaging: Keep in a place protected from rain and direct sunlight.

### **9- INSTALLATION OR CONNECTION TO OTHER PRODUCTS**

The IMPLANT, to be implanted, needs to be connected to a specific key, after being integrated it receives the connection of prosthetic components so that the patient's oral restoration can be carried out.

Attention: The use of any part, accessory or material not specified or provided for in these Instructions for use is the sole responsibility of the user.

#### **10- WARNINGS AND/OR PRECAUTIONS**

For the safe and effective use of Implacil Osstem implants, we strongly suggest that specific training be undertaken, since the surgical techniques required to place these dental implants are highly specialized and complex procedures. Using an improper technique may lead to implant failure and/or loss of supporting bone or other adverse effects.

A rigorous screening of implant candidates should be performed, including age and the presence of parafunctional habits .

Visual inspection and periapical radiographs are essential for determining anatomical accidents, occlusion conditions, periodontal status and bone adequacy.

cephalometric radiography , tomographies

- Do not sterilize or use a second time!
- Do not use an open or damaged package.
- The critical threshold (temperature/time relationship) for bone tissue necrosis is a temperature of 47°C applied for one minute. At a temperature of 50°C for more than one minute, the critical threshold for non-regeneration of peri-implant bone is reached .
- The use of cutting drills, with progressive diameter, with appropriate speed (1500 to 2000 RPM/min), associated with abundant irrigation, allows to avoid excessive thermal elevation.

#### **11-PERFORMANCE ON ESSENTIAL REQUIREMENTS FOR SAFETY AND EFFICACY OF THE MEDICAL PRODUCT, AND ANY UNDESIRABLE SIDE EFFECTS**

##### **11.1 - Indication, purpose or use for which the product is intended**

Indication: Patients who require total or partial oral rehabilitation.

Purpose: To perform the function of the root of the missing tooth

##### **11.2 - Undesirable side effects and contraindications**

Implant installation surgery should only be performed by professionals who have adequate technical training, as it is a highly specialized and complex surgical technique.

Excessive bone loss or implant fracture can occur when the implant is loaded beyond its functional capacity. Physiological and anatomical conditions can impair the performance of dental implants. This should be taken into account when placing dental implants in patients with the following conditions: poor bone quality; poor oral hygiene; medical conditions such as blood disorders or uncontrollable hormonal conditions; metabolic or systemic diseases that affect bone or wound healing, such as uncontrolled diabetes and coagulation disorders; chemotherapy or radiotherapy treatments; untreatable

occlusion or joint problems; cases of insufficient bone height or width; tobacco and alcohol abuse.

Potential adverse events associated with the use of dental implants may include: nonintegration; loss of integration; dehiscence requiring bone grafting; perforation of the maxillary sinus, inferior border, lingual surface or inferior alveolar canal; infection presented by: abscess, fistula, suppuration, inflammation or radiolucency ; persistent pain, numbness, paresthesia; excessive bone loss requiring intervention; and implant fracture.

#### **11.3 - PRODUCT SAFETY AND EFFICACY**

The Implant is completely safe, as long as safety rules are followed in all procedures performed. The product was developed with non-toxic medical-grade material, manufactured in ASTM F-67 titanium. When maintained (stored) and preserved as mentioned in item 8 of these Instructions for use, the product will not lose or alter its physical and dimensional characteristics. The product is sterile and suitable for use for up to 5 years from the date of sterilization.

**Magnetic Resonance Imaging (MRI):** Implacil has not evaluated the safety of its dental implants in the MRI environment. The safety of magnetically induced displacement forces and the extent to which radiofrequency (RF) induced heating may affect patient safety have not been evaluated. The safety of Implacil dental implants in the MRI environment is unknown. Performing an MRI on a patient with this device may result in harm to the patient.

Implacil recommends that medical professionals and the manufacturer of the MRI equipment be consulted for compatibility with implants prior to use. Please be advised that Implacil implants are machined in Titanium Gr.4 ASTM F67.

#### **12-PRECAUTIONS IN CASE OF DISPOSAL**

To avoid environmental contamination or misuse of the IMPLANT, when it is no longer in use, it must be discarded in an appropriate place, that is, in a place for potentially contaminated materials.

#### **13-WARRANTY TERMS**

IMPLACIL OSSTEM - MATERIAL ODONTOLÓGICO SA guarantees the owner of this product a guarantee against any material or manufacturing defect that may arise within 365 days, counting from the date of purchase by the first purchaser.

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---

**9 - Implante Collo HE / Implant Collo HE / Implante Collo HE**

33124 - IMPLANTE COLLO HE 4.0X08  
33125 - IMPLANTE COLLO HE 4.0X10  
33126 - IMPLANTE COLLO HE 4.0X11.5  
33127 - IMPLANTE COLLO HE 4.0X13  
33128 - IMPLANTE COLLO HE 4.0X15  
33135 - IMPLANTE COLLO HE 4.0X18

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Note: The references mentioned here do not accompany the product, must be purchased separately and are subject to specific ANVISA Registration.

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EN:

MA= Open Tray; MF= Closed Tray; AR= Anti Rotational

HI Prosthetic Abutments – Compatible with the HI Cylindrical, HI Conical, HI IL Maestro Conical and Veloce HI Implants.

HE Prosthetic Abutments - Compatible with the Cylindrical HE, Conical HE, Conical HE IL Maestro, Collo HE and Veloce HE Implants.

Morse Cone Prosthetic Abutments (CM) – Compatible with the Morse Cone, Due Cone IL Maestro and Veloce CM Implants.

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EN: These Morse Taper Abutments are compatible with all Morse Taper Implant diameters: 3.3, 3.5, 3.75, 4.0, 4.3, 4.5, 4.75, and 5.0

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EN: These Cone Morse taper are compatible with all Cone Morse Implant diameters: 3.3; 3.5; 3.75; 4.0; 4.3; 4.5; 4.75 and 5.0.

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EN: Implacil certifies that the connection system of the Implants that make up this family is safe and that their combination with other products or equipment as instructed in the instructions for use does not alter the expected performance.

The Implants of this family may only be used in interface with products manufactured by Implacil, whose interface codes are listed in the product's Instructions for use.

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EN: Commercial models - Implant compatibility is determined by diameter, as per the tables below:

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EN:

The Morse Taper implant is supplied with a 1 mm Healing Abutment, a component used to replace the cover screw to protect the prosthetic connection during the healing period and promote soft tissue conditioning. This Healing Abutment is registered with ANVISA under a separate process should the dental surgeon need to purchase it separately

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EN:

\* There is a single installation key that fits all models and its name is CM Installation Key.

\*\* There is a single Transfer for all models and its name is Morse Cone Transfer.

\*\*\* All Morse Cone Prosthetic Abutments fit any implant diameter.

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| Implante Cônico HE / HE Conical Implant / Implante cônico HE | Chave Instalação / Installation Key / Instalación de llaves | Cicatrizadores HE / Healing Aids HE / Cicatrices HE | Transfer HE | Pilares Protéticos HE / Prosthetic Abutments HE / Pilares protésicos HE | Parafuso/Tapa Implante Screw/ImplantCover Tornillo/Cinta de implante |
|--------------------------------------------------------------|-------------------------------------------------------------|-----------------------------------------------------|-------------|-------------------------------------------------------------------------|----------------------------------------------------------------------|
| Ø3,50mm                                                      | Ø3,5mm                                                      | Ø3,5mm                                              | Ø3,5mm      | Ø3,5mm                                                                  | Ø3,5mm Rosca M1.8                                                    |
| Ø4,0mm                                                       | Ø4,0mm                                                      | Ø4,0mm                                              | Ø4,0mm      | Ø4,0mm                                                                  | Ø4,0mm Rosca M2.0                                                    |
| Ø5,0mm                                                       | Ø4,0mm                                                      | Ø5,0mm                                              | Ø5,0mm      | Ø5,0mm                                                                  | Ø5,0mm Rosca M2.0                                                    |

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